

System III Tissue Processing Cassettes

Specially adapted surface optimises print quality

CellPath's System III Cassettes come with an attached lid that has an easy to use 'break-off' hinge. The locking mechanism on the cassettes ensure maximum tissue security during processing.

There are three versions available; slotted, biopsy and micro biopsy, all with a precisely engineered writing surface to elevate the print quality on the most widely used cassette printers.


Our System III Thermal cassettes are optimised for printers using thermal print technology, and the System III Inkjet Cassettes optimised for printers that adopt inkjet technology - *please see table below.*

Printer Type	Printer Model
Inkjet Printer	Leica IPC Sakura AutoWrite
Thermal Printer	Sakura SmartWrite Primera Signature DiaPath AP Itineris

The slotted cassette has a large open surface area to increase reagent flow and improve processing. The biopsy and micro biopsy cassettes incorporate square apertures of 0.9mm x 0.9mm and 0.38mm x 0.38mm respectively, which ensures the safe accommodation of smaller biopsy and micro biopsy specimens, whilst also enabling good levels of fluid transfer.



High quality printing of text and data matrix barcodes.

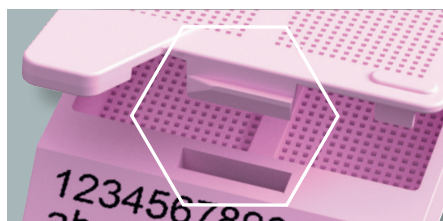

The reprofiled catch improves the opening and closing of the cassette.



LOOSE



THREADED

All System III Cassettes are supplied either loose or threaded for easy loading into the most used cassette printers.



V2 Product UPGRADE



OPTIMISED PRINTING AREA

KEY BENEFITS

- 1 Compatible with the most widely used cassette printers
- 2 Increased reagent flow
- 3 Integral fixed frangible lid
- 4 High rigidity and stability
- 5 Large colour range

Tissue Processing Cassettes

The System III Cassette's New and Improved Features

Improved Specimen Security


- a. Reprofiled catch improves the opening and closing of the cassette, providing an easier and more secure positive closure

Increased Safety

- b. Chamfer frangible hinge minimises any sharp edges when lid is removed

Improved Printing

- c. Increased rigidity during printing ensures a robust and consistent marking surface
- d. Re-engineered surface for optimal clarity of print

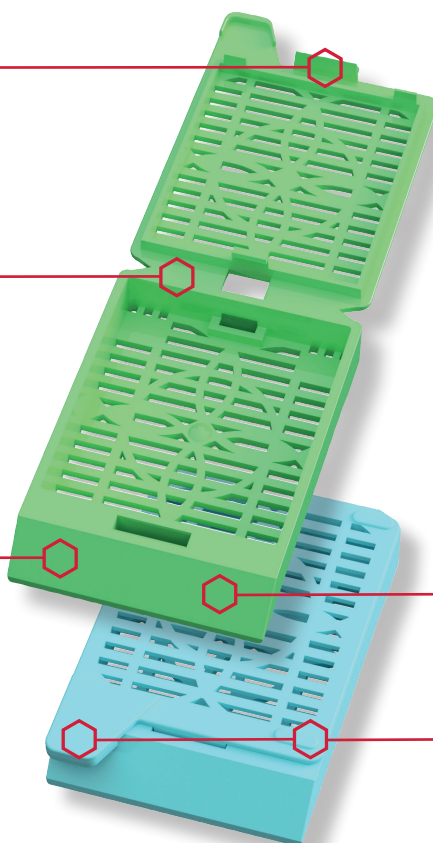

Greater Print Protection


- e. Increased protection against ink being removed during de-waxing on a heated wax trimmer

Increased Stack Stability

- f. 0.25mm raised support points on rear of the lid
- g. Additional 2 upstands to the front of the cassette skirt
- h. Small increase of length to the finger tab

These features help ensure that the cassette is always level and stacked correctly in the printer's hopper. The resultant increase in stability reduces potential downtime in the laboratory.

The System III Cassettes are available in a full range of colours

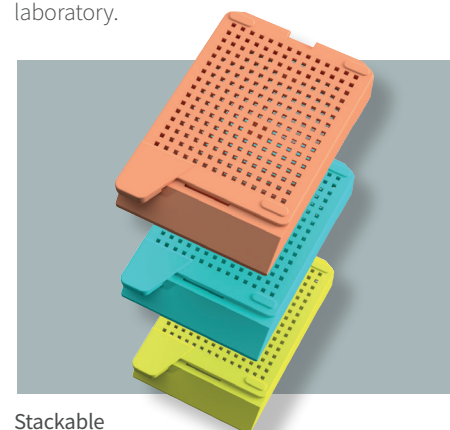

CellPath Ltd

80 Mochdre Enterprise Park,
Newtown, Powys, SY16 4LE, UK

T: +44 (0)1686 611333

E: sales@cellpath.com

cellpath.com



CARING FOR OUR
PLANET

Helping our customers meet their environmental commitments

CellPath are committed to making the world a greener place and keeping our ecological footprint as small as possible. As such, cassette threads are now manufactured from specialist cardboard made from pulp that is recyclable and fully biodegradable.

Tissue Processing Cassettes


Inkjet Printer Compatible System III Slotted Cassettes

Loose		Threaded		Colour
Item Code	Quantity	Item Code	Quantity	
EAI-0202-10A	500	EAI-0202-72B	2000	White 
EAI-0204-10A	500	EAI-0204-72B	2000	Green 
EAI-0205-10A	500	EAI-0205-72B	2000	Grey 
EAI-0206-10A	500	EAI-0206-72B	2000	Blue 
EAI-0207-10A	500	EAI-0207-72B	2000	Orange 
EAI-0208-10A	500	EAI-0208-72B	2000	Pink 
EAI-0209-10A	500	EAI-0209-72B	2000	Yellow 
EAI-0210-10A	500	EAI-0210-72B	2000	Red 
EAI-0211-10A	500	EAI-0211-72B	2000	Peach 
EAI-0212-10A	500	EAI-0212-72B	2000	Tan 
EAI-0215-10A	500	EAI-0215-72B	2000	Brown 
EAI-0216-10A	500	EAI-0216-72B	2000	Dark Blue 
EAI-0218-10A	500	EAI-0218-72B	2000	Dark Green 
EAI-0219-10A	500	EAI-0219-72B	2000	Lavender 
EAI-0220-10A	500	EAI-0220-72B	2000	Turquoise 




Inkjet Printer Compatible System III Biopsy Cassettes

Loose		Threaded		Colour
Item Code	Quantity	Item Code	Quantity	
EAI-0302-10A	500	EAI-0302-72B	2000	White 
EAI-0304-10A	500	EAI-0304-72B	2000	Green 
EAI-0305-10A	500	EAI-0305-72B	2000	Grey 
EAI-0306-10A	500	EAI-0306-72B	2000	Blue 
EAI-0307-10A	500	EAI-0307-72B	2000	Orange 
EAI-0308-10A	500	EAI-0308-72B	2000	Pink 
EAI-0309-10A	500	EAI-0309-72B	2000	Yellow 
EAI-0310-10A	500	EAI-0310-72B	2000	Red 
EAI-0311-10A	500	EAI-0311-72B	2000	Peach 
EAI-0312-10A	500	EAI-0312-72B	2000	Tan 
EAI-0315-10A	500	EAI-0315-72B	2000	Brown 
EAI-0316-10A	500	EAI-0316-72B	2000	Dark Blue 
EAI-0318-10A	500	EAI-0318-72B	2000	Dark Green 
EAI-0319-10A	500	EAI-0319-72B	2000	Lavender 
EAI-0320-10A	500	EAI-0320-72B	2000	Turquoise 

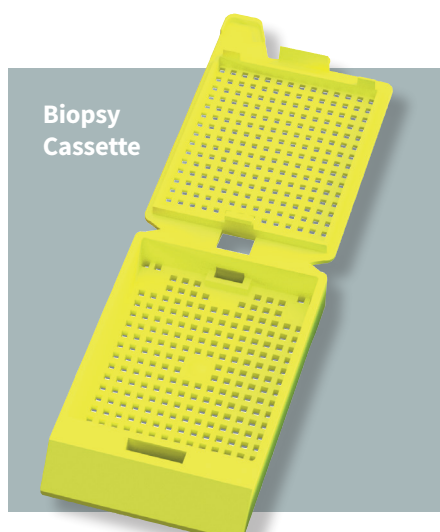

Inkjet Printer Compatible System III Micro Biopsy Cassettes

Loose		Threaded		Colour
Item Code	Quantity	Item Code	Quantity	
EAI-0402-10A	500	EAI-0402-72B	2000	White 
EAI-0404-10A	500	EAI-0404-72B	2000	Green 
EAI-0405-10A	500	EAI-0405-72B	2000	Grey 
EAI-0406-10A	500	EAI-0406-72B	2000	Blue 
EAI-0407-10A	500	EAI-0407-72B	2000	Orange 
EAI-0408-10A	500	EAI-0408-72B	2000	Pink 
EAI-0409-10A	500	EAI-0409-72B	2000	Yellow 
EAI-0410-10A	500	EAI-0410-72B	2000	Red 
EAI-0411-10A	500	EAI-0411-72B	2000	Peach 
EAI-0412-10A	500	EAI-0412-72B	2000	Tan 
EAI-0415-10A	500	EAI-0415-72B	2000	Brown 
EAI-0416-10A	500	EAI-0416-72B	2000	Dark Blue 
EAI-0418-10A	500	EAI-0418-72B	2000	Dark Green 
EAI-0419-10A	500	EAI-0419-72B	2000	Lavender 
EAI-0420-10A	500	EAI-0420-72B	2000	Turquoise 

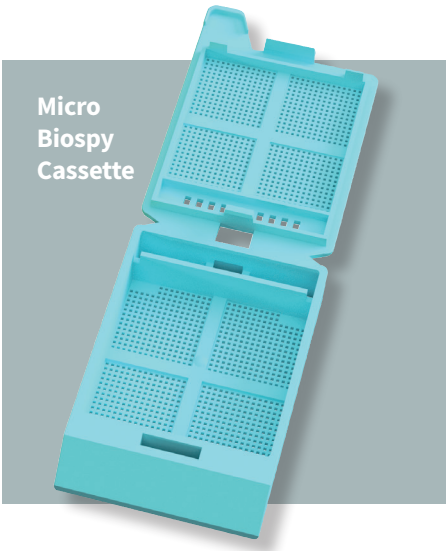

Tissue Processing Cassettes


Thermal Printer Compatible System III Slotted Cassettes

Loose		Threaded		Colour
Item Code	Quantity	Item Code	Quantity	
EAI-0502-10A	500	EAI-0502-72B	2000	White 
EAI-0504-10A	500	EAI-0504-72B	2000	Green 
EAI-0505-10A	500	EAI-0505-72B	2000	Grey 
EAI-0506-10A	500	EAI-0506-72B	2000	Blue 
EAI-0507-10A	500	EAI-0507-72B	2000	Orange 
EAI-0508-10A	500	EAI-0508-72B	2000	Pink 
EAI-0509-10A	500	EAI-0509-72B	2000	Yellow 
EAI-0510-10A	500	EAI-0510-72B	2000	Red 
EAI-0511-10A	500	EAI-0511-72B	2000	Peach 
EAI-0512-10A	500	EAI-0512-72B	2000	Tan 
EAI-0515-10A	500	EAI-0515-72B	2000	Brown 
EAI-0516-10A	500	EAI-0516-72B	2000	Dark Blue 
EAI-0518-10A	500	EAI-0518-72B	2000	Dark Green 
EAI-0519-10A	500	EAI-0519-72B	2000	Lavender 
EAI-0520-10A	500	EAI-0520-72B	2000	Turquoise 



Thermal Printer Compatible System III Biopsy Cassettes

Loose		Threaded		Colour
Item Code	Quantity	Item Code	Quantity	
EAI-0602-10A	500	EAI-0602-72B	2000	White 
EAI-0604-10A	500	EAI-0604-72B	2000	Green 
EAI-0605-10A	500	EAI-0605-72B	2000	Grey 
EAI-0606-10A	500	EAI-0606-72B	2000	Blue 
EAI-0607-10A	500	EAI-0607-72B	2000	Orange 
EAI-0608-10A	500	EAI-0608-72B	2000	Pink 
EAI-0609-10A	500	EAI-0609-72B	2000	Yellow 
EAI-0610-10A	500	EAI-0610-72B	2000	Red 
EAI-0611-10A	500	EAI-0611-72B	2000	Peach 
EAI-0612-10A	500	EAI-0612-72B	2000	Tan 
EAI-0615-10A	500	EAI-0615-72B	2000	Brown 
EAI-0616-10A	500	EAI-0616-72B	2000	Dark Blue 
EAI-0618-10A	500	EAI-0618-72B	2000	Dark Green 
EAI-0619-10A	500	EAI-0619-72B	2000	Lavender 
EAI-0620-10A	500	EAI-0620-72B	2000	Turquoise 

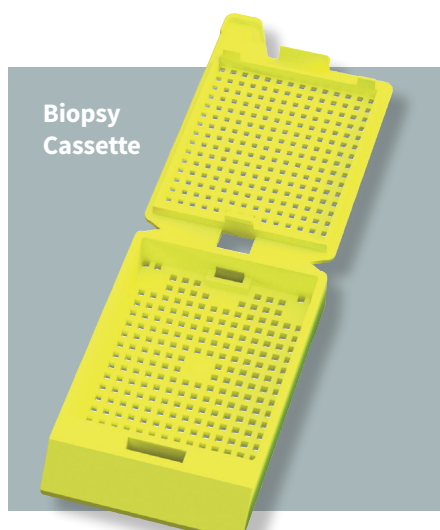

Thermal Printer Compatible
System III Micro Biopsy Cassettes

Loose		Threaded		Colour
Item Code	Quantity	Item Code	Quantity	
EAI-0702-10A	500	EAI-0702-72B	2000	White 
EAI-0704-10A	500	EAI-0704-72B	2000	Green 
EAI-0705-10A	500	EAI-0705-72B	2000	Grey 
EAI-0706-10A	500	EAI-0706-72B	2000	Blue 
EAI-0707-10A	500	EAI-0707-72B	2000	Orange 
EAI-0708-10A	500	EAI-0708-72B	2000	Pink 
EAI-0709-10A	500	EAI-0709-72B	2000	Yellow 
EAI-0710-10A	500	EAI-0710-72B	2000	Red 
EAI-0711-10A	500	EAI-0711-72B	2000	Peach 
EAI-0712-10A	500	EAI-0712-72B	2000	Tan 
EAI-0715-10A	500	EAI-0715-72B	2000	Brown 
EAI-0716-10A	500	EAI-0716-72B	2000	Dark Blue 
EAI-0718-10A	500	EAI-0718-72B	2000	Dark Green 
EAI-0719-10A	500	EAI-0719-72B	2000	Lavender 
EAI-0720-10A	500	EAI-0720-72B	2000	Turquoise 

